




# HIVAMAT<sup>®</sup> 200 Evident

## MODALITY-ASSISTED MANUAL THERAPY

- Pain Relief
- Relaxes Muscles
- Promotes Mobility
- Gentle, Non-Invasive Treatment



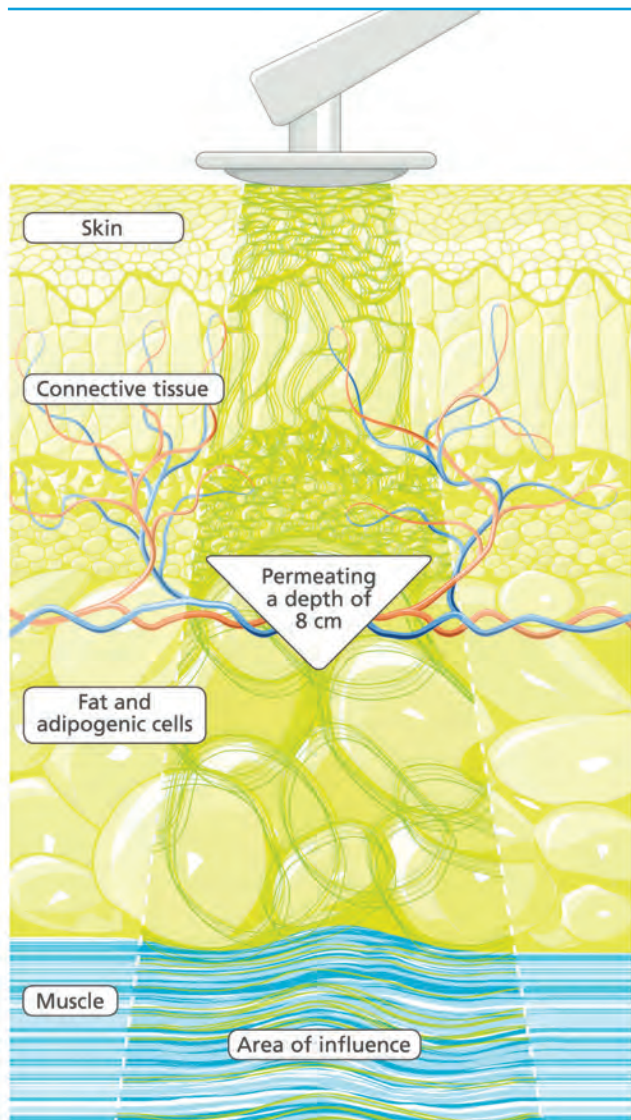

Hivamat 200 Evident is a patented, deep oscillation therapy system. Oscillations of the tissue are created using electrostatic attraction and friction with oscillator heads, resulting in a gentle, yet extremely effective physiological effect on the skin, muscles, connective tissues and lymphatic system. Hivamat allows complete control over treatment parameters with your choice of interface methods: a software-guided digital touchscreen; or via an analog unit with LED display.

## ORDERING DETAILS

02037	<b>A</b> Hivamat® 200 Evident
01004	<b>B</b> Evident Trolley
02036	<b>C</b> Hivamat® 200
379	Applicator handhold
261	Connection cable
778	Mains cable
391	Membrane set 5 cm
392	Membrane set 9.5 cm
600	Operating instructions
381	Oscillator head 5 cm
386	Oscillator head 9.5 cm
260	Patient lead
348	Powder
229	Therapy card
382	Titanium neutral element

## TECHNICAL SPECIFICATIONS

### HIVAMAT 200 Evident

Protection Class	1, Type BF
Power Connection	100 - 240 VAC +/- 10%
Mains Frequency	50 / 60 Hz
Current Consumption	0.2 A at 115 V (0.1 A at 230 V)
Power Consumption	20-30 VA
Output Voltage Max.	400 V
Load Impedance	10 megaohms
Output Frequency	5 - 250 Hz
Dimensions	10.24" x 13.78" x 14.57"
Weight	18.95 lbs



A



B



C